
beurer FC 52 Laguna

Gesichtsbürste | Facial brush | Brosse cosmétique visage | Cepillo facial | Spazzola facciale | Yüz Fırçası

Gerätebeschreibung | Device description | Description de l'appareil | Descripción del aparato | Descrizione dell'apparecchio | Cihaz açıklaması

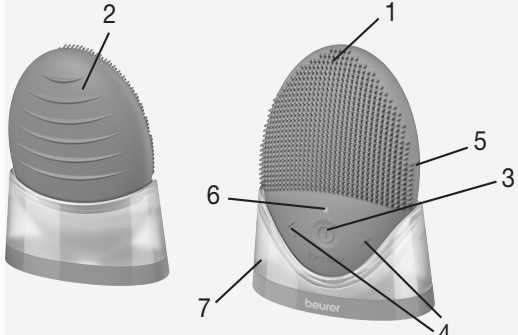- Reinigungsbereich | Cleansing area | Zone de nettoyage | Zona de limpeza | Area di pulizia | Temizleme bölgesi
- Massagebereich | Massage area | Zone de Massage | Zona de masaje | Area di massaggio | Masaj bölgesi
- Ein/Aus-Taste | ON/OFF button | Touche MARCHÉ/ARRÊT | Botón de encendido/apagado | Pulsante ON/OFF | Açma/kapatma tuşu
- Intensitätsregelung - / + | Intensity control - / + | Réglage de l'intensité - / + | Regulación de la intensidad - / + | Regolazione dell'intensità - / + | Yoğunluk kontrolü - / +
- Ladebuchse mit wasserdichtem Verschluss | Charging socket with water-proof seal | Prise de chargement avec fermeture étanche | Toma de carga con cierre impermeable | Presa per ricarica con chiusura a tenuta | Su geçirmez kapaklı şarj konektörü
- LED-Licht zur Akkustands-, Lade- und Intensitätswechselanzeige | LED light for displaying the battery level, charging process and changes in intensity | Lumière LED pour afficher l'état de la batterie, de la charge et du changement d'intensité | Luz LED para indicar el estado de la batería, el nivel de carga y el cambio de intensidad | Luce a LED di indicazione dello stato di carica della batteria ricaricabile, della carica e del cambiamento di intensità | Şarj ve yoğunluk değişimi göstergesi ve pil seviyesi için LED ışık
- Halterung | Holder | Support | Soporte | Supporto | Tutucu

DEUTSCH

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ENGLISH

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Included in delivery

Check that the exterior of the delivery cardboard packaging is intact and make sure that all components are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 x hand-held device
- 1 x USB-C charging cable
- 1 x holder
- 1 x instructions for use

Warnings and safety notes

WARNING

- Danger of suffocation by packaging material. Keep packaging material away from children.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- The hand-held part must be disconnected from the connection cable before it is cleaned in water.
- When using the device in the bathroom, unplug the plug after charging as water in the vicinity constitutes a danger, even when the device is switched off.
- Keep the device away from extended exposure to direct sunlight, extreme heat or boiling water.
- Avoid using corrosive cleaning products or care products which are clay or silicone based.
- The device is only intended for domestic/private use, not for commercial use.
- This device may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely, and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Under no circumstances should you open or repair the device yourself, as fault-less functionality can no longer be guaranteed thereafter. Failure to adhere to this instruction will void the warranty.

WARNING

- Do **NOT** use the device
 - on broken skin,
 - on open wounds,
 - if you suffer from a skin disease or irritated skin,
 - immediately after sunbathing,
 - if you are taking medication that contains steroids.

Intended use

The device is only to be used on the human face (apart from the eye area). The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

Notes on handling batteries

- If your skin or eyes come into contact with fluid from the battery cell, flush out the affected areas with plenty of water and seek medical assistance.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush the rechargeable batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.
- Fully charge the battery prior to initial use.
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

IMPORTANT

Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.

Signs and symbols

The following symbols are used on the device, in these instructions for use, on the packaging and on the type plate for the device:

	Warning Warning instruction indicating a risk of injury or damage to health
	Important Safety note indicating possible damage to the device/accessory
	Product information Note on important information
	Separate the packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	Read the instructions
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.
	Manufacturer
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Do not dispose of batteries containing hazardous substances with household waste
	Direct current The device is suitable for use with direct current only
	Suitable for use in a bath or shower.
	Polarity of d.c. power connector
	Importer symbol

Initial use

Charge the device for at least 2 hours before initial operation. To do this, proceed as follows:

- Insert the mains cable into the charging socket on the side of the facial brush.
- Insert the USB-plug into a suitable mains adapter and connect this with the mains socket.
- Whilst charging, the white LED light above the On/Off button of the facial brush will flash.
- When the battery is fully charged, the LED light will remain on.

Note

If the battery level is very low, the white LED light will light up three times in succession when pressing the On/Off button.

IMPORTANT

Do not use the device while it is charging!

Use

- Moisten the FC 52 facial brush with water. Never use the device while it is dry, as this can irritate the skin.
- Apply the cleansing gel to your facial skin.
- Press the On/Off button once to switch on the device. The LED light will briefly flash once and the facial brush will begin to vibrate.
- The device has a total of 15 intensity levels and when turned on, will begin to vibrate at a moderate speed setting. To switch between the individual intensity levels, press the - / + button. The LED light will briefly flash 1x each time the intensity level is changed.
- Now apply the device gently to your face and cleanse your face using even movements. For cleansing the delicate skin around the eyes, use a lower intensity level. The massage area can also be used for massages.
- After use, switch the device off by again pressing the On/Off button once.
- Cleanse your face thoroughly with water to remove any residue from your skin.
- Pat your face dry with a towel and then apply a moisturising cream using even, circular movements.

IMPORTANT

Stop use immediately if use of the facial brush leads to reactions, and consult your doctor.

Note

To ensure comfortable treatment, do not press the silicone nubs too firmly against the skin. Treat each section of skin for a maximum of 20 seconds.
We recommend a daily cleansing routine of 1 minute mornings and evenings. To achieve an efficient result, duration of use should not exceed 3 minutes. The facial brush switches off automatically after 15 minutes.

Cleaning

IMPORTANT

- Never clean the device with cleaning products containing alcohol, petrol or acetone.
- Clean the device using only the methods specified.
- Thoroughly clean the device after every use. Clean the silicone nubs with soap and water. Afterwards rinse off the soap residue with water.
- Pat dry the device with a lint-free cloth.
- Ensure that no water gets inside the charging socket. If this occurs, only charge the device once the opening has fully dried.

Disposal

- The empty, completely flat rechargeable batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.
- The codes below are printed on rechargeable batteries containing harmful substances:
Pb = Battery contains lead
Cd = Battery contains cadmium
Hg = Battery contains mercury

For environmental reasons, do not dispose of the device in the household waste at the end of its service life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

The device battery must be disposed of separately. The battery must not be disposed of with domestic waste.
To remove the battery, remove the plastic cover at the bottom of the device and then unscrew the housing. Remove the battery from the device and dispose of it at a suitable collection point. For information on disposal, please contact the authorised seller or the public authority of a suitable collection point.

Technical data

Input	5,0 V 1,0 A
Dimensions	11,3 x 7,8 x 4,8 cm
Weight	Approx. 80 g
Battery: Capacity Nominal voltage Type designation	200 mAh 3,7 V Lithium-ion

Warranty

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

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